

Click · Store · Share – Anywhere!

The new generation in video otoscopy



AURICAL OTOcam 300



otometrics

MADSEN · AURICAL · ICS

- Effective and intuitive counseling
- Easy, precise and time-saving diagnostics
- Professional video quality
- Lightweight and portable
- Built for comfort and ease of use
- Powerful light source enables consistent high quality of video and images
- Unique built-in heater
- No rechargeable batteries
- Higher patient satisfaction
- More recommendations





Effective counseling and diagnostics – hand in hand

Personal involvement of the patient is essential for a full understanding of the procedure, whether it is a diagnosis, treatment or hearing instrument selection. Effective counseling increases patient compliance and satisfaction with the solution you recommend. That means more satisfied patients, new recommendations and higher revenues.

The AURICAL OTOcam 300 is a lightweight and portable professional video otoscope for dispensers and ENT physicians. Designed and built for ease of use, comfort and portability, the AURICAL OTOcam 300 video otoscope combines a powerful light source with advanced software to let you display crisp,

full-color video and images on a standard laptop or desktop computer or even a large screen. And the built-in heating system dramatically reduces annoying delays caused by condensation.

Take it even further – let OTOcam 300 be the tool for training your patients in how to operate and maintain their hearing instrument. Showing tiny details on a large screen is just so much more efficient. This training will save your time by reducing unplanned patient support.

- Fast and simple installation
- Easy to learn
- Easy to use
- Allows you to spend your time efficiently
- Aesthetic, space-saving and portable design
- Ergonomic, lightweight handpiece for comfort and confidence
- Specula with and without cerumen management

Now including specula with cerumen management





Plug and play and ease of use gives you high efficiency

The AURICAL OTOcam 300 video otoscope is easy to install, easy to learn and easy to use. You can set up the AURICAL OTOcam 300 video otoscope in a matter of minutes and start using it immediately. That means you can start offering your patients a new level of involvement and service from the moment you unpack and set it up.


Easy-to-handle specula ensure optimum hygiene. Available with or without cerumen management, the specula are easy to snap on and off. The new specula increase efficiency of the patient flow, and support professional service for hearing care professionals. *

* Local restrictions related to carrying out cerumen management services may apply.

The lightweight and ergonomic design of the AURICAL OTOcam 300 takes otoscopy to new places, letting you bring video otoscopy out to your patients – you don't have to bring a separate computer – install the software, plug in and you're ready!

Let OTOcam 300 give you that competitive edge.

- Intuitive software
- Variety of viewing modes and annotation functions
- NOAH-compatible for complete office integration
- Effortless sharing of reports, data, video and images
- Reporting and data exchange
- Integrated control buttons

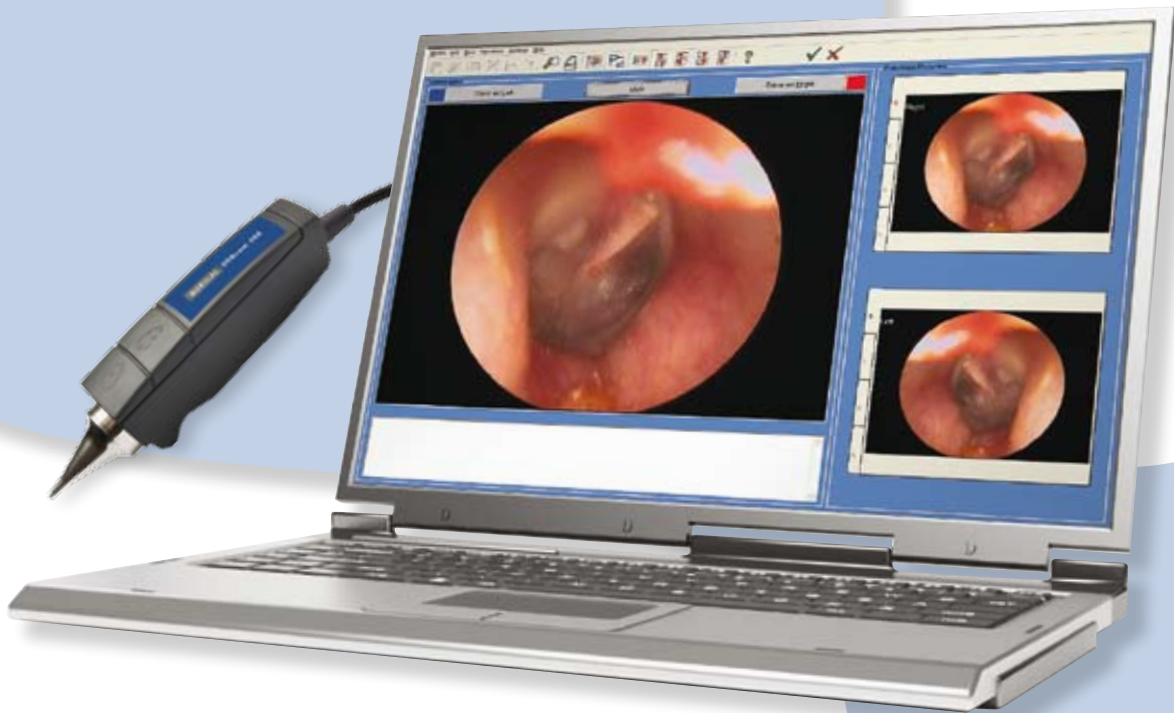



Software integration – for your benefit

The intuitive AURICAL OTOcam 300 software module allows you to do your job faster, smarter and with fewer errors. With or without NOAH, the AURICAL OTOcam 300 software will provide you with annotation and reporting functionality that meets today's requirements for paper-based and electronic reporting.

Unique integration between the software and AURICAL OTOcam 300 allows you to focus on your primary tasks without worrying about the equipment. All you have to do is aim, click the button on the camera, and let the software store the image without even touching the computer!

The flexible compatibility of the AURICAL OTOcam 300 means you can also use it with other standard video editing software to produce illustrative video clips.

Video otoscopy has never been easier or more effective!

For more information or to order your AURICAL OTOcam 300 video otoscope for delivery as soon as possible, contact your local Otometrics representative.

Technical data

Optical data

- Sensor resolution 752 (H) x 582 (V) pixels
- Minimum working distance 8 mm

Preheating system

- Fog-free after 5 min.

Electrical data

- Supply voltage through USB port

Connectors

- 2x USB, type A, USB 2.0

Mechanical data

- Weight including cable 320 g /0.7lbs
- Cable length 2950 mm/9.6 ft.

Safety

- Degree of protection IP 67
- Application class BF
- Protection class II DIN EN 60601-1

Computer requirements

- 1.5 GHz Pentium® IV CPU (or better depending on individual operating system requirements)
- Minimum 1GB RAM and 64 MB graphics memory
- Windows® XP (professional) with SP1 or higher, and Windows Vista 32 bit with SP1
- 2 powered USB 2.0 ports for AURICAL OTOcam 300 video otoscope
- CD-ROM drive
- NOAH 3.1, NOAH 3.1 for ENTs, or higher for NOAH mode operation

Looking for associated products?

MADSEN Itera II



AURICAL Plus



HI-PRO USB



AURICAL Visible Speech



Contact your local
Otometrics representative

GN Otometrics, Europe. +45 45 75 55 55. info@gnotometrics.dk
GN Otometrics, North America. 1-800-289-2150. sales@gnotometrics.com
www.otometrics.com