



Vestibular *precision*

The VF5 Video Frenzel system from Interacoustics provides a cost effective means of recording and observing fixation free eye movements on patients suffering from balance disorders – right in your own office.

Advantages of the VF5

- Lightweight Mask
- Fixation-free observation
- Documentation on VCR/DVD-R
- Observation of rotational eye-movement

The lightweight mask is the most important component of the VF5 Video Frenzel system. Since the mask provides the main patient interface it must be comfortable and non-threatening to the

patient. With this in mind, great care was taken in the design and engineering of the VF5 mask in order to provide a comfortable fit for any patient, even young children. The main advantage of any frenzel mask is that it completely prevents patients from visually fixating and thereby suppressing spontaneous nystagmus which is a critical indicator of vestibular disorders.

The added benefit of the VF5 is that the 100% light occlusion provides a total freedom from visual fixation and the examiner can easily view enlarged images of the patients eye movements on the connected monitor. Even rotational eye-movements that are normally difficult to observe and record during a routine

exam are evident with this frenzel system.

Another benefit of the VF5 is that you can record your observations by simply attaching a VCR or DVD-R. Making a patient recording then allows for a convenient pre and post treatment comparison.

Monocular or Binocular

The VF5 can be purchased with either one infra-red camera (monocular) or two infra-red cameras (binocular). The binocular system provides the convenience of immediate side by side comparisons of the left and right eye movements on your patients and a reduction in evaluation time.

VF5m



Camera




Combi Mask



Camera Control Box

Not included

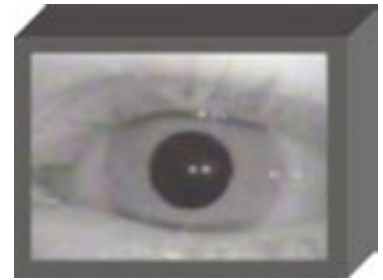

VCR/DVD-R and/or a variety of standard TV screens can be connected

VF5b



Camera



Combi Mask




Camera



Camera Control Box

Not included

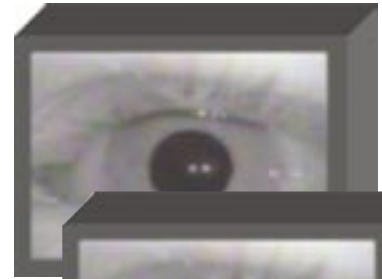

VCR/DVD-R and/or a variety of standard TV screens (one or two) can be connected

83002101-US - 3 - 08/2006

Interacoustics USA

Phone: 1-800 947 6334

Fax: 1-952-903-4100

E-mail: info@interacoustics-us.com

Web: www.interacoustics.com

Mail: 7625 Golden Triangle Drive

Eden Prairie, MN 55344

